## MFP NANO PLUS 4000





The MFP Nano plus 4000 is specially designed to clearly determine the separation efficiency of HEPA and ULPA filter media in accordance with DIN EN 1822-3 and ISO 29463-3. The U-SMPS is a modern and powerful nanoparticle measuring device with a measuring range of 5 nm to 1  $\mu$ m for particle size and quantity analysis. In the MFP Nano plus 4000, the separation efficiency at a specific size is measured simultaneously with one UF-CPC condensation core counter each in raw and clean gas.

#### 工作原理

# SIMULTANEOUS DETERMINATION OF THE FRACTIONAL SEPARATION EFFICIENCY WITH TWO UF-CPCs

The simultaneous measurement of the fractional efficiency in the MFP Nano plus 4000 offers the following special advantages:

- The real-simultaneous measurement of the particle concentration in raw and clean gas halves the measurement time when determining the fractional efficiency.
- The combination of two UF-CPC versions, the UF-CPC 200 for highest concentrations up to 2,000,000 particles/cm³ (single counting mode) in the raw gas and the UF-CPC 50 for highest counting rates at low concentrations in the clean gas corresponds to a dilution factor of 1:200. Conventional aerosol dilution can therefore be dispensed with for highly efficient filter media.


The universally applicable aerosol generator UGF 2000 can be used to generate aerosol distributions with DEHS or salt (NaCl / KCl) that are tailored to the MMPS range.

The extensive automation of the test sequence together with the clearly defined individual components and the individually adjustable sequence programs of the FTControl filter testing software ensure the high reliability of the measurement results.

The MFP filter test stand is a modular filter test system for flat filter media and small mini filters. Pressure loss curve, fractional separation efficiency or loading can be determined reliably and therefore economically within a very short time.

MFP Nano plus 4000



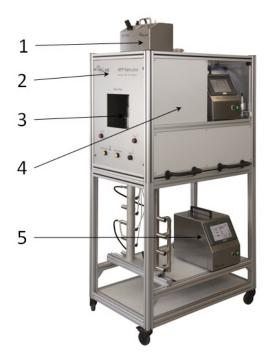

- 1. Variable aerosol generation with the UGF 2000 for KCl / NaCl or DEHS. Integrated Nafion drying section: The volume flow for aerosol generation is controlled individually via mass flow controller.
- 2. Aerosol neutralization: Soft X-ray Charger XRC 049 or corona discharge (optional)
- 3. Mobile pneumatic filter holder: This is used to quickly remove and load the test stand with integrated DEMC 2000 (Differential Electrical Mobility Classifier). The DEMC 2000 (Differential Electrical Mobility Classifier) contained in the U-SMPS pre-classifies the polydisperse aerosol generated with the UGF 2000 in terms of particle size. After the DEMC 2000, only monodisperse particles are contained in the test aerosol. The corresponding size is automatically regulated by the DEMC 2000 control unit.
- 4. Raw gas measurement of the number distribution with UF-CPC 200 for high concentrations: To measure the number of particles in the raw gas, the condensation core counter UF-CPC 200 is used for particle counting up to 2,000,000 particles/cm<sup>3</sup> in single counting mode. This means that no dilution system is required for raw gas measurement in very high concentrations. There is no need for time-consuming cleaning of the systems.
- 5. Clean gas measurement of the number distribution with UF-CPC 50 for low concentrations: The "full flow" UF-CPC 50 is optimized for low particle concentrations (single counting mode up to 10,000 particles/cm<sup>3</sup>) in clean gas. With "full flow" particle counting in the UF-CPC, the entire sampling volume flow is fully analyzed. This means that very high counting rates are achieved in the clean gas at low particle concentrations.

The Palas filter testing software FTControl controls the U-SMPS and evaluates the data.

By selecting the appropriate solution concentration, the particle size distribution generated in the MFP Nano plus 4000 is adapted to the respective MPPS range.

版本: 2025年10月16日 MFP Nano plus 4000



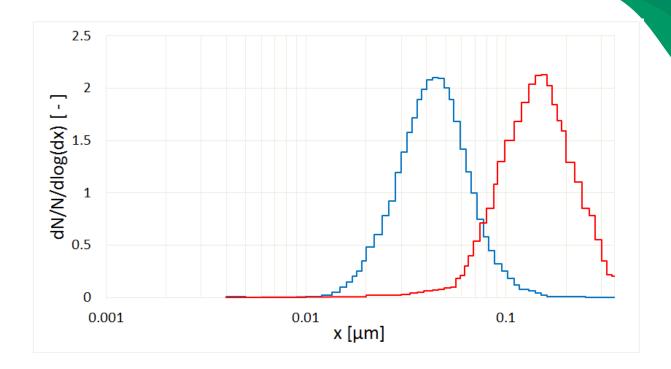


Abb. 1: Adaptation of the particle sizes to the required MPPS range with DEHS

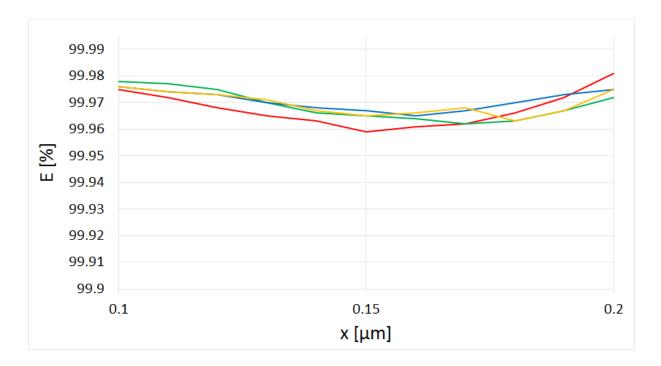


Abb. 2: Comparison of fractional separation efficiencies in the MPPS range at 140 nm

#### Automation:

The MFP Nano plus 4000 has built-in mass flow controllers for volume flow control, which can be automatically controlled via the FTControl filter testing software. The sensor data such as volume flow and differential pressure at the filter are automatically recorded during the filter test.

版本: 2025年10月16日 MFP Nano plus 4000



Proof of conformity of sampling in raw and clean gas:

Proof of compliance of the raw and clean gas measurement without a filter is demonstrated during the briefing.

With the MFP Nano plus 4000, fractional efficiency measurements in the MPPS range as well as over the entire measuring range are easily possible. In addition, the respective pressure loss of the medium is clearly determined at the respective flow velocity.

版本: 2025年10月16日 MFP Nano plus 4000



### 优势

- Real-time determination of the fraction separation efficiency above 20 nm
- Half the measurement time thanks to simultaneous measurement in the upward and downward directions
- No dilution necessary!
- $\bullet$  Combining two UF-CPC versions, the UF-CPC for the highest concentrations, up to 2,000,000 particles/cm<sup>3</sup> (single count mode) in the raw gas and the UF-CPC 50 for top counting rates at low concentrations
- Simple use of different test aerosols, such as NaCl / KCl or DEHS (others on request)
- · Simple measurement of the fraction separation efficiency and determination of the MPPS range
- · High reproducibility of the test method
- Cleaning can be performed independently by the customer
- Short set-up times, fast throughput times
- Mobile set-up, easy to move on castors

### 标准和证书

EN 1822-3, ISO 29463-3



## 技术数据

气溶胶	Dusts (e.g., SAE dusts), salts (e.g., NaCl, KCl), liquid aerosols (e.g., DEHS)
滤材测试面积	100 cm <sup>2</sup>
测量范围(粒径)	U-SMPS: 10 – 800 nm
体积流量	$0.48 - 5.76  m^3/h$ - pressurized operation
电源	115 – 230 V, 50/60 Hz
Differential pressure measurement	0 – 2,500 Pa (others on request)
Inflow velocity	1.3 – 16 cm/s (others on request)
Compressed air supply	6 – 8 bar
Dimensions	Approx. 760 • 2,100 • 985 mm (H • W • D)

MFP Nano plus 4000



## 应用领域

- Testing of filter media and small mini filters in product development and production monitoring
- Testing capability according to DIN EN 1822-3 (HEPA / ULPA) and ISO 29463-3
- Fractional efficiency measurement for other filter media in the range of approx. 20 nm to 1  $\mu$ m



Mehr Informationen:

https://www.palas.de/zh/product/mfpnanoplus4000